
# MANUAL DE INSTRUCCIONES

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### CA-01 LIGHT

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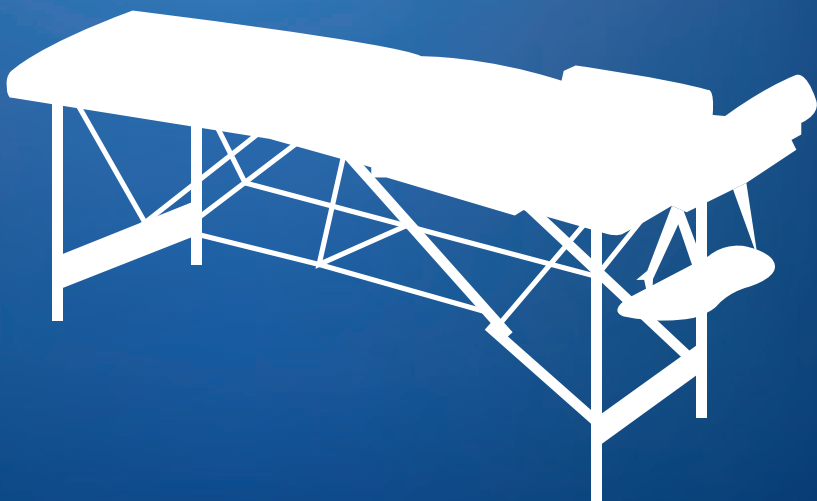## IMPORTANT

Keep the instructions for future reference. Read these instructions before assembling and using this product. Adult assembly is required.

Thank you for choosing this aluminium massage table CA-01 LIGHT! In this instruction manual you will find all the information you need to assemble, use and maintain this product.

Before you start, we recommend that you read this instruction manual carefully to ensure safe and optimal use of the product.

## 1. PRODUCT CHARACTERISTICS

<b>NAME</b>	Aluminium massage table
<b>MODEL</b>	CA-01 LIGHT
<b>REFERENCE</b>	M-00061/08-CR

## 2. ASSEMBLY INSTRUCTIONS

1. First loosen the buckle. Place the stretcher on its side and unfold it.
2. Open the legs. Leave the lower legs in their original position when you start to unfold the stretcher.
3. Straighten both parts of the stretcher.

4. Lift the stretcher until all four legs touch the floor.

5. Open the head restraint. Insert the bars of the head restraint into the holes and fit the pillow.
6. Loosen the side handle, adjust the height and angle of the head restraint. Then push the handle back.
7. Attach the armrests and check both sides of the stretcher before use.

### 3. STRETCHER STORAGE

1. Place the stretcher on its side and fold down the legs completely.
2. Close both halves of the stretcher by pressing and locking the locking catch.

### 4. WARNINGS

- The stretcher should be folded with the help of a person.
- Check and tighten all screws before use.
- Always place the product on a flat and stable surface. To prevent the stretcher from destabilising, make sure that all legs are at the same level, otherwise injuries may occur.
- Make sure that all screws are tightened securely and that the cables are not worn or damaged. Also check that there are no cracks anywhere on the stretcher.
- Do not overload the stretcher, it must not exceed the recommended weight limit as this could cause serious injury to the user. Do not place fingers in the recess of the stretcher.
- Do not use alcohol-based products to clean the stretcher cover as this may cause damage to the material. Use a damp cloth only.

### 5. WARRANTY CONDITIONS

This product is covered by a 36 month warranty against manufacturing defects.

Thank you for choosing Mobiclinic! We appreciate your confidence in our products. We are committed to providing you with the best quality and service.

We know you will enjoy your new CA-01 LIGHT as much as we enjoyed creating it for you. If you have any questions or need assistance, please do not hesitate to contact us at the email address below: **clientes@mobiclinic.com**.

We are here to help you.

**ES** - Capacidad máxima del producto 250 kg  
**EN** - Maximum product capacity 551,15 lbs  
**IT** - Capacità massima del prodotto 250 kg  
**FR** - Capacité maximale du produit 250 kg  
**PT** - Capacidade máxima do produto 250 kg  
**DE** - Maximale Produktkapazität 250 kg  
**NL** - Maximale productcapaciteit 250 kg  
**SWE** - Maximal produktkapacitet 250 kg  
**PL** - Maksymalna pojemność produktu 250 kg  
**DK** - Maksimal produktkapacitet 250 kg

Made in P.R.C.

HANGZHOU FIRSTWELL INDUSTRY CO., LTD.  
Room 708, Building 2, Zheshang Caifu Center, 97#  
Gudun Road, Hangzhou, Zhejiang, China  
anny@hzrstwell.com

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GRUPO R. QUERALTÓ S.A.