Manual de instrucciones

TO-01

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- | | |
|------------------------------------|--|
| 1. Introduction and classification | 7. Settings |
| 2. Product features | 8. Measuring method |
| 3. Technical structure | 9. Replacement of batteries |
| 4. Shape structure | 10. Maintenance |
| 5. Display icon definitions | 11. Information of electromagnetic compatibility |
| 6. The function keys | |

1. INTRODUCTION AND CLASSIFICATION

Notice: Notice to the user and/or patient: any serious incident that has occurred in relation to the device must be reported to the manufacturer and to the competent authority of the Member State in which the user and/or patient is established.

Indications: The infrared thermometer is suitable for persons of all ages.

The TO-01 device allows body temperature to be taken by receiving infrared energy radiating from the forehead. It is hygienic, accurate and easy to use. The user simply directs the infrared probe to the forehead, presses the measure button and the body temperature is displayed on the LCD screen.

The device is commonly used in hospitals and as an aid to daily living.

It can also be used to measure temperatures from 0°C to 118°C on objects.

The device should not be used in oxygen, air or flammable gas, or nitrogen oxide combinations.

2. PRODUCTS FEATURES

- High precision infrared sensors
- Adaptability to ambient temperature. Can be used in complex environments.
- Patented sensor designed to ensure more accurate measurement.
- Fever alarm function.
- Ability to store 32 temperature readings.
- Measurement units: Celsius (°C) or Fahrenheit (°F).
- Automatic shut-off and power saving functions.

3. TECHNICAL PARAMETERS

| | | |
|--------------------------------|---|---------------------------------------|
| Temperature measurement method | | Non - contact |
| Effective distance | | 5cm - 15cm |
| Measurement Mode | Body mode | 32.0°C - 42.0°C (85.9°F - 109.2°F) |
| | Surface mode | 0°C - 118°C (32 °F - 244°F) |
| Measurement Accuracy | Body mode | ±0.2°C / 0.4°F |
| | Surface mode | ±1.0°C / 1.8°F |
| Display resolution | 0.1°C / °F | |
| Operating environment | 16°C - 35°C / 61°F - 95°F (≤80% moisture condensation) | |

OBJECT OF SURFACE

| | |
|-----------------------|--|
| Storage condition | -25°C - 55°C / -13°C - 131°F (≤90% moisture condensation) |
| Electrical source | DC3V (2 AA Batteries) |
| Power consumption | When off 10mw |
| | When measuring ≤30mw |
| Power level indicator | Indication for low power level |
| Remember have | Before measuring the 32 sets of data |
| According to the unit | Celsius (°C) o Fahrenheit (°F) |
| Auto power-off | In 10 seconds |
| Dimensions | 88 mm * 44mm * 153 mm (length x weight x height) |
| Weight | 130 g (including battery) |

4. SHAPE STRUCTURE

1. Infrared sensor
2. LCD
3. Backlight On/Off button
4. Up button
5. Down button
6. Function key and unit of measurement conversion
7. Power on
8. Battery compartment cover
9. Handle

5. DISPLAY ICONS DEFINITIONS

| Icon to define | Icon | State description |
|----------------------------|---|--|
| Sound indication |  | Sound prompt (buzzer) |
| | | Buzzer to shut down, no prompt |
| Measuring mode | Body temp | Body temperature mode |
| | Surface temp | The object temperature measurement model |
| Temperature unit | °C | Celsius degree |
| | °F | Fahrenheit degree |
| According to the numerical | 888.8 | Temperature value |
| Stand in one's own light |  | Backlight LCD display |
| Mnemonic symbol | LOG | Record the latest 32 sets memories |

ICONS

| | |
|---|--------------------------------------|
|  | Recycling garbage |
|  | Warning |
|  | Flammables and explosives |
|  | Read the operation manual before use |
|  | Class BF medical device |

6. THE FUNCTION KEYS

| Function key | Functional description |
|----------------|---|
| Backlist key | Increase the backlight function, for the use of the night |
| Option key | Adjust the settings of parameters |
| Conversion key | Hold down to select temperature units °C/°F into settings |

7. SETTINGS

The thermometer is provided with 4 types of settings:

- Temperature units.
- Sound switch.
- Temperature deviation.
- Temperature taking mode: The mode setting is carried out on the mode change button.

The other settings must be made from the main menu of the thermometer.

The menu is as follows:

| Menú | Functions | "-" | "+" | Initial value | Remarks |
|------|-------------------------|---------------------------------|-------------------------------|---------------|-----------------------|
| F1 | Temperature alarm point | Reduce to 0.1°C | Increase of 0.1°C | 37.8°C | Scope of 37°C-42.9°C |
| F2 | OffsetTC | The downward migration of 0.1°C | The upward migration of 0.1°C | | The range of -2°C-2°C |
| F3 | Indications switch | Close | Open | Open | Press +/- to select |
| F4 | Measuring mode | Surface Object | Human body | Human body | Press +/- to select |

TEMPERATURE ALARM SETTINGS: F1

Press the "MODE" key, the display will show "F1". From the set alarm, press the "+" button to increase the temperature by 0.1°C or the "-" button to decrease the temperature by 0.1°C. To finish, press the "MODE" button to confirm.

TEMPERATURE COMPENSATION SETTINGS: F2

Press the "MODE" key, the display will show "F1", press MODE again to display "F2". Then set the temperature compensation. Finally, press the MODE key again to save the update you have made.

WARNING SETTINGS: F3

Press MODE key, the display will show "F1", press twice more to display "F3", press "MODE" again to set the tone, the display will show "ON". Press "+" or "->" to make the display show "OFF" and disable the prompts, finally press "MODE" again to save the update you have made.


MEASUREMENT MODE: F4

Press "MODE" key, the display will show "F1", press "MODE" key three times to display "F4", press "MODE" again to choose the temperature taking mode, then the display will show "surf", press "+" or "-" and the display will show "body". To save the selected temperature mode press "MODE".

8. MEASURING METHOD

TAKING BODY TEMPERATURE

- Press the measurement key to turn on the device, then the display will show the backlight and all the symbols will appear on the display. The device will beep twice when the self-test is finished and the display will show the last value saved in the memory and then enter the measuring status.
- Make sure that the measurement mode is set to body temperature mode.
- Align the thermometer probe with the forehead, the distance should be about 3-5cm. Press the measurement key. The temperature will then appear on the display when a signal sounds and 0.5 seconds elapse. If the temperature exceeds the limit (the default value is 37.8°C), an alarm will sound.
- When the measurement is finished, OFF will appear on the display and it will turn off after 30 seconds.

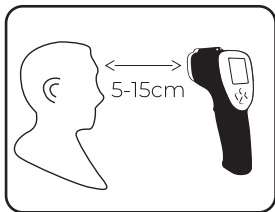NOTE:

1. Keep the sensor and internal cavity clean before taking the temperature.
2. Use the thermometer in a stable temperature environment. If the ambient temperature changes greatly, wait for about 10 minutes before taking the temperature.
3. Do not measure body temperature after measuring extremely high or low temperatures, wait 10 minutes before measuring.

TAKING THE TEMPERATURE OF THE OBJECT

- Press the measurement key to turn on the device and activate the measurement status.
- Make sure that the upper display is not set to body temperature mode, but to surface temperature mode.
- Align the thermometer sensor with the centre of the object, the distance should be 5 to 15 cm. Press the measurement key, the temperature will be shown on the display after a beep, in 0.5 seconds.
- The thermometer will automatically switch off after 10 seconds inoperative.

9. REPLACEMENT OF BATTERIES

- Battery life is approximately 3000 readings in one year.
- When the thermometer indicates low battery, it is recommended to replace the batteries.
- Open the cover and remove the old batteries.
- Insert the new batteries and be careful with the position of each battery.

TIPS

- When not using the thermometer for a while (about 3 months), please remove the batteries.

10. MAINTENANCE

- Mantenga los sensores y la cavidad interna limpios. De lo contrario, las mediciones podrían dar resultados inexactos.
- This device requires only routine cleaning, at intervals that depend on the user's needs.
- To clean the surface, use a clean, soft cloth or cotton swab dipped in sanitary alcohol or water to remove dirt.
- To clean the sensor and internal cavity, use a clean, soft cloth or cotton swab with a little medical alcohol to gently wipe the sensor and internal cavity. Do not use the thermometer until it is completely dry.
- The process of use does not require any maintenance.

11. TROUBLE SHOOTING

| The diagnosis | Disposal measures |
|------------------------------|---|
| The screen shows "Lo" o "Hi" | <ol style="list-style-type: none">1. Check the object on which the temperature is to be taken; note that the accuracy of the temperature measurement on the forehead is not guaranteed if the device is placed on hair, water, perspiration or cosmetics.2. Check the temperature deviation setting.3. In case the temperature change is too large or in case you use the thermometer with an object that is at low temperature after using it at high temperature, it is recommended to use it after 10 minutes to achieve a good thermal equilibrium.4. Check the recommended distance, which should be 5cm to 15cm. |

| | |
|-------------------------------------|--|
| Buttons have no response | <ol style="list-style-type: none"> 1. To load and upload the battery. 2. Check wether the operation under the settings menu. |
| No display or display abnormal | To load and upload battery. |
| Check weather voice to set up close | Check the battery, loading and unloading battery again |

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SHENZHEN BRAV ELECTRONIC TECHNOLOGIES CO., LTD
MedNet GmbH
Grupo R. Queralto S.A.